



LAPAROSCOPIC FORCEPS



VECTEC® SINGLE USE LAPAROSCOPIC DEVICES




MONOPOLAR SINGLE USE LAPAROSCOPIC GRASPING FORCEPS



REF : * MD5 330



REF : F5 330

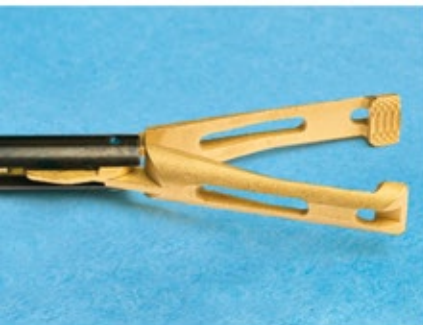

REF : ECS5 330

Reference	Type of jaws	Working length	Diameter	Jaws size
* MD5 330	● Dissector Maryland	330 mm	5 mm	20 mm
A5 330	● Atraumatic	330 mm	5 mm	18 mm
ECS5 330	● Endoclinch	330 mm	5 mm	22 mm
F5 330	● Fenestrated	330 mm	5 mm	18 mm
G5 330	● Claw, 2 x 6 teeth	330 mm	5 mm	18 mm

* without ratchet ● Yellow: Dissector ● Green: Atraumatic ● Purple: Traumatic

Ergonomic handle with ratchet mechanism guarantees strong and precise control with minimal effort during grasping and holding tissues.

Double-action jaws allows perfect grip.

Jaws rotates 360° and are adjusted using the rotational knob at the front of the handle.

The grasping forceps comes fitted with a standard monopolar connector pin.

Electrical insulation is total right up to the active portion of the jaws.

Packaging : box of 5 individually packaged units.

MD Class IIa and IIb
Sterilized by ETO
See instructions for use



Manufactured by:
Péters Surgical
42, rue Benoît Frachon
93013 Bobigny Cedex - FRANCE
Tel: +33 (0)1 48 10 62 62
Fax: +33 (0)1 48 91 22 99
commercial@peters-surgical.com
www.peters-surgical.com